reesa®

Sonic toothbrush

TSA8015P

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CS

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DE

Owner's manual

EN

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Thank you for purchasing TEESA appliance. Read this instruction manual and keep it for future reference. Producer does not take responsibility for damages caused by inappropriate handling and use of the product.


SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Before connecting the device to power supply socket, make sure the voltage indicated on the device corresponds to the voltage in the power supply socket.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon sound when using
 - before cleaning
- When you disconnect the power cord of the device, grasp and pull the plug, not the cord/cable.
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Fixed power cable may only be replaced in an authorized service point.
- Do not use this device if it malfunctions or has been water damaged. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.
- Keep the device and its cord/cable away from heat, water, moisture, sharp edges and any other factor which may damage the appliance or its cord.
- Do not use the product for purposes different from displayed in this operating manual.
- Producer of this product is not responsible for damages caused by inappropriate handling and using of the device.
- Disconnect from power outlet if the device is not going to be used for a long time.
- Use only authorized accessories.
- Do not use the device in the tub or shower.
- Do not use this device in the place, where it can fall into water.
- In such case, disconnect the power cord from the power adapter.
- Do not handle this device with wet hands.
- Appliance for household use only. This device is not intended for professional use.
- NOTICE: The device has waterproof handle, and is secured from moisture (can be used in the bathroom), but do not submerge it underwater.

PRODUCT DESCRIPTION

1. Brush head
2. Handle
3. Power button
4. Change mode/intensity button
5. Intensity indicator
6. Clean mode indicator
7. Soft mode indicator
8. Strong mode indicator
9. White mode indicator
10. Polish mode indicator
11. Power/charging indicator
12. Charger
13. USB power cord

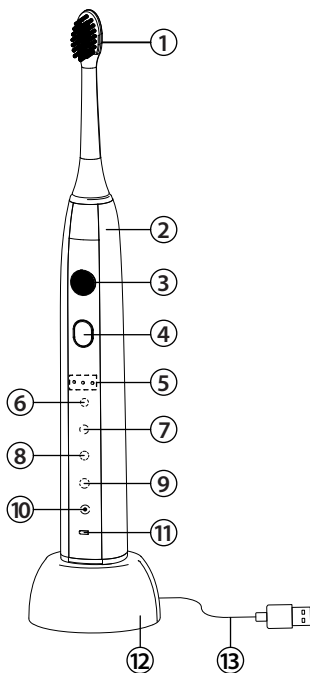CHARGING

1. Before first use, make sure, that the battery of the device is fully charged.
 2. Connect the USB plug to the power adapter.
 3. Place the toothbrush in the charger.
 4. Power/charging indicator will glow red, indicating charging in progress.
 5. Full charging cycle takes about 16 hours and allows for up to 30 days of usage.
 6. Once the charging is completed, the power/charging indicator will glow green, indicating full battery level.
 7. Disconnect the charger after finished charging.
- Note:** if the battery level is low, the power/charging indicator will flash red.

OPERATION

1. Before using, attach a brush to the handle and apply toothpaste.

Note: if a brush has been installed correctly, there should be 0,5 mm of space between the brush head and the handle.

1. Before turning on the toothbrush, select the work mode by pressing the change mode/intensity button.
2. To avoid splashing, guide the brush-head to your teeth before switching the appliance on.
3. Place the bristles against the teeth at a slight angle towards the gum line.
4. Press the power button to turn it on.
5. During brushing, there are three intensity levels to choose from. To change intensity level, press the change mode/intensity button.

TIMER

1. The device has been equipped with auto switch off function:

- 2 minutes in modes: Clean, Soft, Strong, White.
- 1 minute in Polish mode.

2. The toothbrush will stop shortly every 30 seconds to remind user the passing time.

Tip: This feature allows to measure equal time for each tooth quadrant. This way you'll spend 30 seconds on each part.

EN

MODES

This device offers different brushing modes:

1. Clean - removes tartar and bacteria from teeth, high frequency
2. Soft - for sensitive teeth and gums, gentle frequency
3. Strong - removes tartar build up, high frequency
4. White - whitens teeth, clean and polish modes switching alternately
5. Polish - polishes and brightens teeth, high strength and frequency

- The device enters Clean mode when powered on for the first time. Later on, the toothbrush will turn on in the recently used mode.
- Selected work mode indicator will glow on the toothbrush.

BRUSHING

- Guide the brush slowly and gently from tooth to tooth. Do not scrub teeth as you would with traditional toothbrush.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth. It is suggested to brush the teeth in order: outside of the teeth, insides, chewing surfaces.
- Do not press hard or scrub, let the brush do all the work.

CLEANING

- After brushing, rinse the brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Wipe them dry before completely before you put the handle onto the charging unit.
- Before cleaning the charger, make sure to unplug it from power mains. Clean the charger only with soft, slightly damp cloth.

IMPORTANT: Do not clean the charger by submerging it in the water or under running water.

- After cleaning, place the charger on a flat surface.
- Do not plug the charger to the power adapter if it is still wet.

SPECIFICATION

- 5 work modes: Clean / Soft / Strong / White / Polish
- 3 intensity modes
- Strokes per minute in work modes:
 - Clean: 32000
 - Soft: 31000
 - Strong: 36000
 - White: 31000-36000
 - Polish: 46000
- IP rating: IPX7
- Auto power-off after 2 minutes, 30 second timer
- Auto power-off after 1 minute, 30 second timer (only in Polish mode)
- Medium bristles
- Noise level: \leq 68 dB
- Power / charge indicator
- Inductive charger
- Usage on one charge: up to 30 days
- Battery: 800 mAh
- Power supply: 110-240 V; 50/60 Hz
- In set: 2x brush head, charging unit, travel case

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for LECHPOL ELECTRONICS Sp. z o.o. Sp.k., ul. Garwolińska 1, 08-400 Miętne.

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