teesa®

Oral irrigator

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Owner's manual EN

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SAFETY INSTRUCTIONS

Thank you for purchasing TEESA appliance. Read this instruction manual and keep it for future reference. Producer does not take responsibility for damages caused by inappropriate handling and use of the product.


In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Before connecting the device to power supply socket, make sure the voltage indicated on the device corresponds to the voltage in the power supply socket.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before cleaning
- When you disconnect the power cord of the device, grasp and pull the plug, not the cord/cable.
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Fixed power cable may only be replaced in an authorized service point.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.
- Keep the device and its cord/cable away from heat, water, moisture, sharp edges and any other factor which may damage the appliance or its cord.
- Do not use the product for purposes different from displayed in this operating manual.
- Producer of this product is not responsible for damages caused by inappropriate handling and using of the device.
- Disconnect from power outlet if the device is not going to be used for a long time.
- Use only authorized accessories.
- Do not use the device in the tub or shower.
- Do not use this device in the place, where it can fall into the water.
- Appliance for household use only. This device is not intended for professional use.
- NOTICE: The device has waterproof handle, and is secured from moisture (can be used in the bathroom), but do not submerge it underwater.

- Attention: when you take the device out of the box for the first time, there might be some liquid on the bottom of the water tank. It is distilled water, used in the factory to test the water flow of the device. It is not dangerous and is not a warranty claim. Pour out this water and rinse the water tank in the sink.

CHARGING


1. Before first use, make sure, that the battery of the device is fully charged.
2. Connect the USB plug to the power adapter.
3. Connect the charging plug to the irrigator.
4. The mode / charging indicator will glow red, while it is charging.
5. Full charging cycle takes about 4 hours and allows for up to 45 minutes of usage.
6. After finished charging, the charging indicator will glow blue, indicating that the battery is fully charged.
7. Disconnect the charger after finished charging.

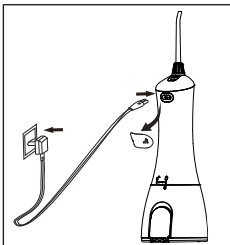
Attention: if the battery level is low, the mode / charging indicator flashes red.

OPERATION


1. Open the water inlet and fill the water tank with regular water. Water temperature cannot be higher than 40°C. Close the water inlet.

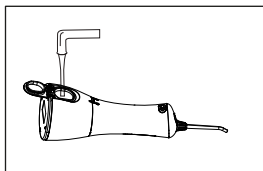2. Before first use, pour a small amount of water to the nozzle slot.

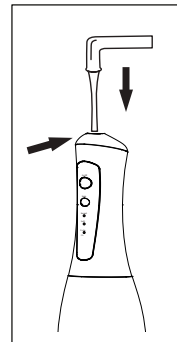
3. Put the nozzle to the irrigator, until you hear a click.

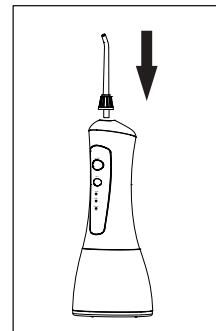4. Press the Mode button to select the work mode.
 - Normal: for everyday irrigation.
 - Soft: for sensitive gums.
 - Pulse: for cleaning and massaging gums.
5. Put the nozzle into your mouth and point it towards the teeth (make sure to keep a small distance between nozzle and teeth, do not put the nozzle directly in front of the teeth).
6. Slightly open your mouth, so the water may pour out.
7. Press the ON/OFF button to turn on the irrigator.


- a. Point the water stream at a 90° angle to your teeth and gums.
 - b. Move the water stream along the teeth line.
8. The irrigator will automatically turn off after 2 minutes. To turn off the irrigator before that time, press the ON/OFF button.

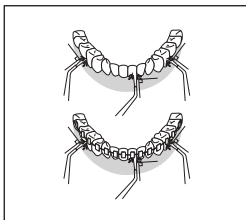After using


1. After finished work, press the nozzle release button and take out the nozzle from the irrigator.
2. Twist the water tank to open position (open padlock symbol) and detach the water tank from the irrigator. Pour out the water from the water tank.
3. Press the ON/OFF button to turn on the irrigator and pump the rest of the water out. Turn off the irrigator.
4. Wipe the irrigator with a piece of cloth.

Attention: if the irrigator is not going to be used for a long time, wipe it carefully and make sure that it is dry all the way.

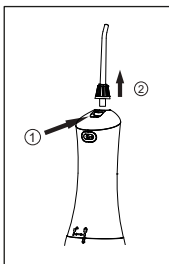CLEANING

- Clean the bigger stains with soap and water.
- To clean the device, do not submerge it in water.
- After cleaning the device, dry it completely.
- Nozzles must be replaced every 6 months.

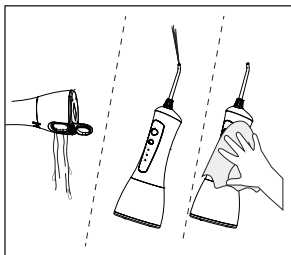SPECIFICATION

- Precise cleaning of gums and between teeth
- 2 minute timer
- 3 work modes:
 - Normal (6,2 bar)
 - Soft (3,1 bar)
 - Pulse (2,7 - 6,2 bar)
- 3 tips included: standard, delicate, orthodontic
- Power: 5 W
- Water flow: 180 ml/min
- Protection class: IPX7
- Noise level: <75 dB
- Water tank capacity: 278 ml
- Working time: approx. 45 minutes
- Charging time: approx. 4 hours
- Battery: 3,7 V / Li-ion / 1500 mAh
- Charger: 5 V / 1 A
- Weight: 244 g
- Dimensions: 79 x 210 x 89 mm (without tip)
- In set: irrigator, charger, USB cable, tips

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this material from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for LECHPOL ELECTRONICS Sp. z o.o. Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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