BEAUTIFLY

Glossy Acne

Young
series

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
Technical data:

Size:	161x161x44.5mm
Weight:	110g
Rated Power:	5W
Rated current:	1A
Battery:	3.7V polymer lithium ion battery
Battery capacity:	500mAh
Included:	satin bag, USB cable, user manual

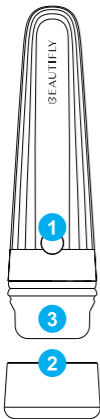**Safety statement**

Read the instructions carefully before use. Follow these safety guidelines when using the device.

Part description

1. Power button
2. Protection cap
3. Treatment surface

Warnings and cautions

- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not immerse the device in water.
- When plugged in, do not leave the device unattended.
- Do not use the device if the electric cord or power supply has been damaged or is working improperly.
- Use the device only with the electric cord and charger that come with it.
- Keep the device out of reach of children.
- Make sure your hands are completely dry when operating the device.
- Turn the device OFF after use or before cleaning and maintenance.
- Follow the instructions for using skin care products during treatment.

Safety instructions:


1. Use this device only for its intended purpose, as described in this manual.
2. Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas. Do not use the device in and around the eye and lip area.
3. If you experience any discomfort while using the device, discontinue use immediately.

How to use

Charging

1. Make sure that the battery is **fully** charged before first use.
2. When charging, insert the **electrical cord** into the charger **plug** and then connect it to the general household power supply. To charge the product, the other end of the electrical cord is inserted into the charging hole in the device. During charging, the LED indicator emits **light**, which goes off when the battery is **fully** charged.

The charging time of the device is about 1.5 hours and its working time is about 2.5 hours.

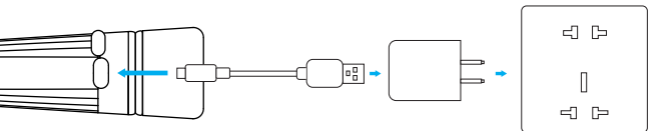Power on/Power off

Power ON: Short press the POWER button to turn on the device and start the **Blue LED Photon Light** mode.

Power OFF: Short press the POWER button to turn off the device.

Low battery reminder

When the product battery is **low**, the **red light** at the button **will flash** once every 2 seconds.

Timing function

The timing function is 3 minutes (the **buzzer will prompt** once every 5 seconds).

Method of use

1. **Clean** the skin **thoroughly** before treatment.
2. Short press the power button to turn on the device.
3. Point the device at the acne area and scan it for five seconds. Move the device over the entire treatment area, making sure that **all** parts affected by acne are scanned **sequentially**.
4. Short press the power button to turn off the device.
5. Take an appropriate amount of skin care cream and **apply** it **evenly** to the acne areas. **Gently** massage the product over the face.
6. **Clean** the treatment surface of the device with a damp cloth or paper towel and cover it with the protection cap.

Note

It is recommended to perform a skin test 24 hours before the first use to ensure that there are no adverse reactions. Do not use the device on and around the eye area.

Recommended treatment guidelines

- It takes 3 minutes to treat the entire face. For shorter treatment time, focus on specific areas of concern. The treatment time for each person is different.
- To achieve optimal results, use the device 1-2 times a day for a period of 1-2 weeks.
- To achieve long-term results after finishing the full course of treatments, we recommend using the device 2-4 times a month or as needed.

Warnings

- Do not use the device if you have a pacemaker or internal defibrillator.
- Do not use the device if you have metal implants in the treatment area, such as gold threads. However, this does not apply to dental fillings or implants, which pose no problem.
- Do not use the device if you are currently diagnosed with cancer, have a history of cancer, especially skin cancer or premalignant moles.
- Do not use the device if you suffer from severe conditions, such as cardiac disorders.
- Pregnant or breastfeeding women should not use this device.
- Do not use the device if you are suffering from immunosuppressive diseases such as AIDS, or are using immunosuppressive medication.
- Do not use the device if you have any disease stimulated by heat, such as herpes simplex recurring in the treatment area.
- Do not use the device if you have any active skin condition within the treatment area, such as, but not limited to, sores, psoriasis, or a rash.
- Do not use the device if you have a history of keloid scarring, abnormal wound healing or fragile skin.
- Do not use the device if you have had a surgical procedure, laser resurfacing, or deep chemical peeling in the treatment area within the past three months, or are still undergoing a process to complete healing.
- Do not use the device if you have permanent fillers within the treatment area.
- If you have temporary dermal fillers within the treatment area, consult your physician before using the device.
- Do not use the device if you have had botox administered in the treatment area within the past three months.
- Do not use the device if you have suspicious lesions within the treatment area.
- Do not use the device if you are currently taking medications, herbal preparations, vitamins, or food supplements that might make your skin dry or fragile.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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EN: The manufacturer declares on their own responsibility that the Glossy Acne Young product complies with the European Parliament and Council Directive 2014/30/EU of February 26, 2014, on the harmonization of member states' legislation concerning electromagnetic compatibility, as well as with the European Parliament and Council Directive 2011/65/EU of June 8, 2011, on the restriction of the use of certain hazardous substances in electrical and electronic equipment, along with Commission Delegated Directive (EU) 2015/863 of March 31, 2015, amending Annex II to European Parliament and Council Directive 2011/65/EU with regard to the list of restricted substances. The full text of the declaration of conformity is available at the following website address:

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