
BEAUTIFLY

OptiRelax Pro Smooth

FitRecovery
series

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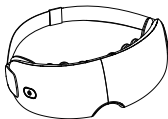
 /beautifully.eu  @beautifully.eu  beautifully.eu

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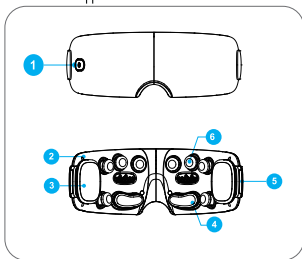- www.beautifly.eu
- www.youtube.com/beautiflyeu

Rated voltage:	5 May
Rated power:	5 W
Rated frequency:	50/60 Hz.

Safety instructions:

1. The product should be stored in a dry and shaded place, at room temperature, out of the reach of children.
2. Never use your device if it is damaged or has a damaged charger, cord or plug.
3. Never charge if the product and charger are immersed in water.
4. The product must not be used in case of eye diseases (after eye surgery, glaucoma, cataract, detached retina and other serious eye diseases).
5. The product is prohibited for use by people who have pacemakers.
6. Do not disassemble the internal components of the product by yourself to prevent the product from malfunctioning.
7. It is prohibited to clean the product with water or corrosive liquid.
8. Before turning on, apply the product to your eyes. If you wear contact lenses, hair ornaments or other items, you should remove them.
9. The product does not work while charging.

Device description [page 2-f1]:

1. Enable/disable.
2. Foundation.
3. Pneumatic massage head.
4. Metal massage head (cooling/heating mode).
5. Support strap.
6. Point vibrating massage heads.

Troubleshooting:

1. The product does not work: Check whether you have pressed the on/off button.
Check if the battery is not exhausted. Check if the time is set to default (press on/off, charge device, restart).
2. Cannot charge: Check if the plug is inserted correctly.
3. The indicator light is dim or flashing, the air pressure is too weak: charge the device.

A description of the action:

1. The default working time after power on is 15 minutes,

the power will be turned off automatically after this time, if you want to continue it, just press the OFF/ON button again.

2. Press and hold the OFF/ON button to turn on the product.

3. Double-click the OFF/ON button to enable silent mode. Double-click again to exit the mode.

4. OFF/ON key: short press - switch modes in order:
Mode 1: Eye heating mode (vibration + air pressure + heat)

Mode 2: Ice Mode (Vibration + Air Pressure + Cooling Mode)

Mode 3: Vitality mode (vibration + air pressure + alternating heat and cold)

Bluetooth:

1. Find EYE-MASSAGER in Bluetooth function and click to connect.

2. After Bluetooth connection has been established once, turn Bluetooth on again and the connection will be established automatically.

3. Power off: Press and hold the switch key for 2 seconds, the power button light will turn off.

4. Blue light - the product is not connected to Bluetooth.
Purple light - Bluetooth connection is on.

Landing:

1. The icons glow red when the product is charging. Fully charged device - the green light is on.

2. It takes about 3 hours to fully charge the battery, and the light is always on white when fully charged.

3. Low battery reminder: when the power is less than 20%, the light is red.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under his sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of legislation. Member States in the field of electromagnetic compatibility and Directive 2011/65/EC of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Commission Delegated Directive (EU) 2015/863 of 2011. On 31 March 2015, amends Annex II to Directive 2011/65/EU of the European Parliament and of the Council on the list of restricted substances.

Model: OptiRelax Pro Smooth.