
BEAUTIFLY

Głowica - MicroGlow
0.5mm

P R O

series

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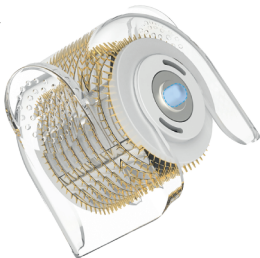Safety Information:

1. The needle cartridge must be moved continuously and can be held on the same skin area for 3-5 seconds, no more than twice.
2. After using the device, a temporary redness may appear on the skin. If these symptoms do not disappear over time, consult a doctor.
3. Avoid skin exposure to the sun, use of spicy foods, seafood, alcoholic beverages and other stimulating foods for three days after treatment.
4. The needle cartridges can only be used once and must be safely disposed of after use.
5. The needle cartridges are for personal use only. Do not share them with others to avoid infection.
6. Do not use the needle cartridge if the disposable packaging is damaged or cracked.
7. Perform a skin test on your hand before first use. Do not use the device if your skin is sensitive or has an allergic reaction.
8. Parts of the device other than the needle tips must not come into contact with liquid to avoid damaging the device.
9. The needle cartridge must be used by one person only.
10. In the event of a failure, do not attempt to repair or disassemble the device yourself.
11. Do not use on the eyelids, thyroid area, lips, ears, breasts or genital area.

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautify Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Głowica - MicroGlow 0.5mm