# BEAUTIFLY

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

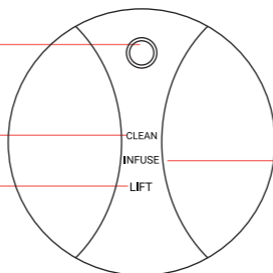
RF connector

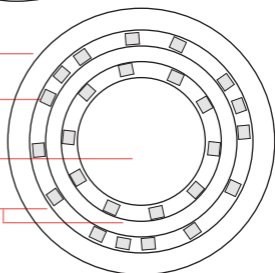
on/off button  
mode switch

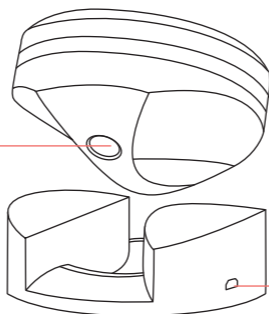
Put the beauty device into the charging base. Three battery charging indicators will flash red when device is charging. After charged, the charging indicator turn green for 5 minutes, then all lights turn off.

(Note: If the red lights not flash, please put the B-Flash PRO in the opposite direction.)

## ► Accessories

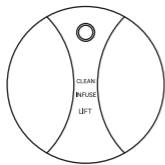B-Flash PRO

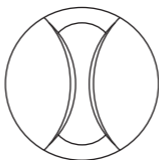wireless  
charging base

USB charging  
cable

## ► Parameters

size	70x40mm
weight	129g
material	ABS plastik, stain


## ► Function

1. Ultrasounds is a high-level sound wave technology that penetrates deeper than the dermal layer. In the skin it increases the temperature of the surrounding tissue, thus increasing oxygenation to the cells, stimulating the regeneration of collagen protein and aiding in the removal of waste products.
2. Green light (520nm wavelength): brighten complexion, sooth and calm skin, reduce oil secretion and balance the water and oil.
3. Red light (630nm wavelength): bio-active light can open pores to help better absorption, improve the active cells, brighten the skin tone, and resist aging and oxidation.
4. Blue light (415nm wavelength): it has the improve effect of removing acne, repairing skin, leaving no scars and pigment
5. Ion-: Negative ions with assistance of high frequency vibration can effectively promote the absorption of nutrition essence. improve skin status from inside to outside, and achieve deep moisturizing effect
6. EMS (Electrical Muscle Stimulation): Through the periodic change of the micro-current to consume fal, compact skin and shape V face contour.
7. RF (Radio Frequency): The main action is overheating collagen fibers to a temperature of about 54 ° C. Under the influence of heating, they shrink and, secondarily, also tension. Overheating contributes to the stimulation of fibroblasts, and thus the recovery and production of new collagen and elastin.  
Thanks to this phenomenon, the effect of increased skin tension, improvement of its elasticity, density as well as improvement of oval and face color is obtained. The B-Flash PRO device uses a frequency of 0.5 MHz

## ► 1. Clean


Before starting the treatment, make-up facial cleanser should be done and the make-up remover should be thoroughly cleaned. We recommend using the Beautifully sonic brush for this purpose. Only properly prepared skin is able to absorb nutrients from cosmetic preparations.
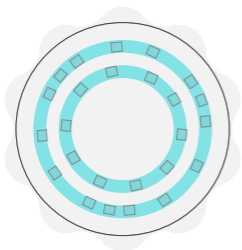
### Step 1

Ultrasonic waves combined with the positive ion to promote lymphatic drainage, and stimulate collagen production, while blue LED light technology kills bacteria for a deep facial detox.

**Note:** You can fix the makeup removal cottons with the stainless steel ring to remove makeup with makeup removers.
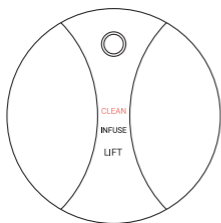
### Step 2


Switch to **CLEAN** mode, apply moisturizing gel primer or serum, and glide across entire face. For optimal results, repeat for 5 minutes, particularly focusing on any areas with redness or acne.

### Step 3


Place the device in parallel to the cheek (as a starting point), then gently rotate from the outside to the chin, from the bottom up to the forehead. Repeat 3 times.

## ▶ 2. Infuse
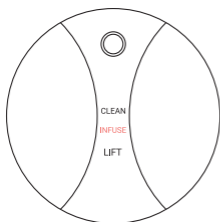
### Step 1

Ion infusion, combined with low-frequency waves and sonic pulsations, enhances absorption of active ingredients in products, leaving your skin supple, while green LED light targets redness and hyperpigmentation.
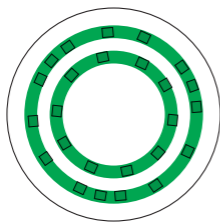
### Step 2



Switch to **INFUSE** mode, apply serum, and glide across face. As product is pushed deep into dermal layers, you will be left with an enviable glow. Use until product is fully absorbed into skin.

### Step 3


Apply serums evenly on the face and place the device in parallel to the chin (as a starting point), gently lift the face from the inside to outside.

### ▶ 3. Lift
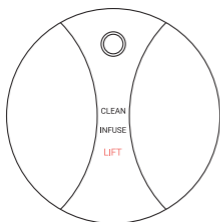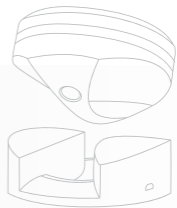
#### Step 1

Micro-currents, RF, and red LED light technology activate your skin's cellular rejuvenation process, increase collagen production, and reduce fine lines and wrinkles.
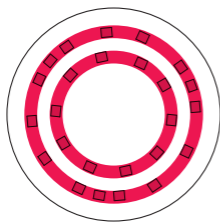
#### Step 2

Switch to **LIFT** mode, apply moisturizing gel primer, and glide across face for 5 minutes, focusing on areas with fine lines and wrinkles.

**\*\*A conductive gel or serum must be used with this step\*\***

#### Step 3

Apply conductive gel evenly on the face and rotate the device from the inside to outside, from the bottom up to the forehead.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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## TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address [serwis@profit-plus.eu](mailto:serwis@profit-plus.eu)  
The exact warranty conditions are under the address:

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I represent that I have read and accepted the terms and conditions specified in the Warranty Certificate.

---