# REMINGTON®

Compact

# CONTROL

..... HPL

**IPL3500**

[www.remington-europe.com](http://www.remington-europe.com)

1

## HAIR COLOUR


## SKIN COLOUR


2

3

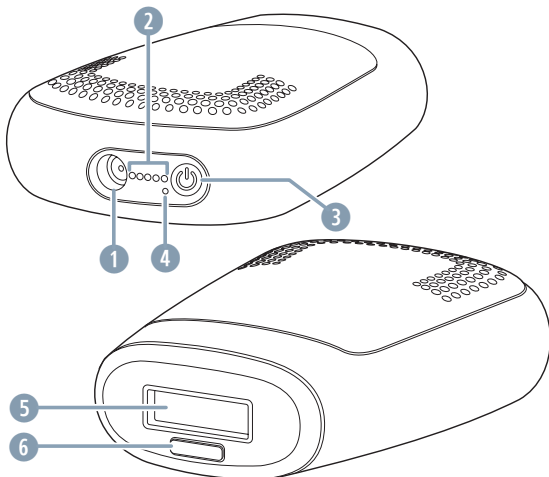4

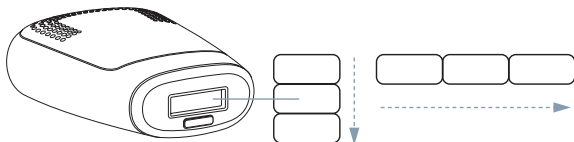The Compact Control HPL is a light-based device for permanent hair reduction designed for home use. It is a powerful electrical device and should be used with special attention to safety. Be sure to read all warnings, contraindications and safety precautions before use, and follow them strictly. We recommend you re-familiarise yourself with the instructions before each treatment.

### INTENDED USE

The Remington Compact Control HPL is intended to remove unwanted hair and effect long-term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs regrowing after a treatment regime. You may use the Remington Compact Control HPL to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face (female facial use only), back, shoulders and chest.

---

### IMPORTANT SAFEGUARDS

---

- This appliance should not be used by anyone under the age of 18 and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge unless they have been supervised/instructed and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by persons under 18.
- Keep the appliance and cable out of reach of children under 8 years.
- Do not use the product with a damaged cord. A replacement can be obtained via our International Service Centres.
- ⚡ The body of this appliance is not washable or water resistant. Do not put the appliance in liquid; do not use it near water in a bath tub, basin or other vessel; and don't use it outdoors.
- Possible hazardous optical radiation emitted from this product. Do not stare at the light source which may be harmful to the eye.
  - Use this appliance only for its intended use as described in this manual.
  - Do not use the appliance if it is damaged or malfunctions.
  - Do not leave the appliance unattended while plugged in.
  - Keep the power plug and cord away from heated surfaces.
  - Only use the parts supplied with the appliance.

- Do not twist or kink the cable, or wrap it round the appliance.
- Do not plug or unplug the appliance with wet hands.
- This appliance is not intended for commercial or salon use.

## **SAFETY WITH YOUR DEVICE**

Home Pulsed Light (HPL) technology in your device provides extra safety with lower energy levels. HPL achieves long-term hair reduction results at a fraction of the energy level used in other light-based hair reduction equipment. The HPL technology reduces the risk of harm or complications, and contributes to your overall safety. Your device has additional features to protect your eyes and your skin.

- **Eyes:** Light pulses are only active when the treatment window is in full contact with the skin. Nevertheless, do not look at the treatment window when you are using the device. This may cause damage to the eyes.
- **Skin:** Your skin is never in direct contact with the optical lens through which light pulses are emitted.

## **BEFORE YOU START USING THE COMPACT CONTROL HPL**

**Be sure to read all warnings and safety information. Before you begin, check to see if the device is suitable for you.**

Use the skin and hair chart at the front of this manual and the integrated skin sensor to determine if this device is right for you (Fig 1).

**WARNING:** The Compact Control HPL is not suitable for very dark or black skin. It is not effective on light or grey hair and is less effective on blonde and red hair.

For facial use: Only suitable for female facial hair on the cheeks, sideburns, chin, neck areas and above the lips (Fig 2).

DO NOT use on the eyes, eyebrows or eyelashes. Doing so can cause serious and permanent eye injury.

## **SKIN TYPE**

- Do not use on naturally dark skin, as it may result in burns, blisters and changes in skin colour.

## **EXPOSURE TO THE SUN**

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for 4 weeks before and 2 weeks after your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.
- To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen, SPF 30 or higher, for 2 weeks after each hair reduction session.

## **AREAS NOT TO TREAT:**

- Do not use on the scalp or ears.
- Do not use on face or neck (male).
- Do not use above the cheekbone (females).

- Do not use on nipples, areola or genitals.
- Do not use if you have tattoos or permanent makeup in the area to be treated.
- Do not use on dark brown or black spots such as moles, birth marks or freckles.
- Do not use on an area of recent surgery, deep peel, laser resurfacing, scars or skin that has been damaged with burns or scalds.

#### **WHEN NOT TO USE/WHEN TO AVOID USING COMPACT CONTROL HPL**

- DO NOT use if you have tanned skin in the intended treatment area, as your skin may be too dark and treatment in tanned areas could temporarily cause skin swelling, blisters, or could change the colour of your skin. Wait until your tan has faded before treating.
- DO NOT use on skin with cosmetics applied, as treating in those areas could temporarily cause skin swelling, blisters or could change the colour of your skin. Make sure to thoroughly clean the intended treatment area prior to use.
- DO NOT flash the same treatment spot multiple times during the same session, as this could temporarily cause skin swelling, blisters, or could change the colour of your skin.
- DO NOT treat any area more than every two weeks, as this will not speed up results, but could increase likelihood of skin swelling, blistering or changing the colour of your skin.
- DO NOT use on damaged skin, including open wounds, cold sores, sunburns, or areas where you have had a skin peel or other skin resurfacing procedure within the last 6-8 weeks.
- DO NOT use on irritated skin, including rashes or swollen skin. These conditions can make the skin more sensitive. Wait for the affected area to heal before use.
- DO NOT use if your skin is sensitive to light. In order to determine if you can use the device, perform a test on a small patch of skin and wait 24 hours to determine if there are any adverse reactions before treating more broadly.
- DO NOT use if you are or may be pregnant or are breast-feeding. This device has not been tested on pregnant women. Hormonal changes are common during pregnancy or breast-feeding and could increase skin sensitivity and the risk of skin injury.
- DO NOT flash against any material other than skin.
- DO NOT use this device for any purpose other than hair reduction.
- DO NOT treat any area where you may someday want your hair to grow back. The results may be permanent and irreversible.
- DO NOT use on the face if you have a history of cold sores; use of HPL or any treatment that irritates skin can cause a recurrence of cold sores.
- If you've had professional or clinical or self-administered laser or IPL or electrolysis hair removal treatments, wait at least 6 months before using this device. If you've waxed, plucked or epilated hair in the treatment area, wait at least 6 weeks before using. The device uses the hair under your skin to disable your hair at the root, so it is not effective if the hair has been pulled out.
- Do not use over or near anything artificial such as silicon implants, Implanon contraceptive implants, pacemakers, subcutaneous injection ports (insulin dispenser) or piercings.

- DO NOT use directly over or near Active Implantable Medical Devices such as implantable pacemaker, internal defibrillator, or insulin devices. Although no problems associated with the use of this device have been reported during studies, it is recommended to observe all active implantable medical device warnings regarding the use of electronic devices to prevent any failure or interruption of an active implantable device.
- Wireless communications equipment may affect the device and should be kept at least 3.3 m away.
- Do not use over or near fillings or dental implants.

## PARTS (FIG 3)

- |                                      |                        |
|--------------------------------------|------------------------|
| 1. Adapter socket                    | 5. Treatment window    |
| 2. Energy level indicator lights (5) | 6. Skin Tone Tester    |
| 3. Control button                    | 7. Adaptor (not shown) |
| 4. Ready light                       |                        |

## INTRODUCTION

The efficiency of hair reduction varies from person to person according to body area, hair colour, and how the device is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair reduction.

A typical hair reduction plan during a full hair growth cycle includes:

- Treatments 1-6: Plan two weeks apart.
- Treatments 7-8: Plan four weeks apart.
- Treatments 9+: Treat as needed, until desired results are achieved.

## EXPECTED TREATMENT OUTCOMES

Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

Skin Type <sup>1</sup>	Body Hair <sup>2</sup>	Area	Average Treatment Time	Effectiveness of Hair Reduction
I - V	Dark blonde to brown	Lower/upper leg	Session 10 - 12 approx. 4 weeks apart	60%
		Arms		50%
		Bikini line		60%
		Underarms		55%

Skin Type <sup>1</sup>	Body Hair <sup>2</sup>	Area	Average Treatment Time	Effectiveness of Hair Reduction
I - V	Brown to black	Lower/upper leg	Session 8 - 10 approx. 4 weeks apart	70%
		Arms		60%
		Bikini line		70%
		Underarms		65%

<sup>1</sup>Skin types are based on the Fitzpatrick scale of skin colours. Refer to the skin colour chart at the front of this manual (Fig 1).

<sup>2</sup> Not effective on white body hair.

### ENERGY LEVELS

The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair reduction results, as well as the risk of possible side effects and complications.

- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 energy level indicator lights.
- When the device is turned on, the energy level is set to level 1 automatically. Only one energy level indicator light will be on.
- To change the energy level, press the control button repeatedly to increase the energy level from 1 to 5. Once at level 5, press the control button again to reduce the energy level from 5 back to 1.
- To determine the right energy level for your treatment area, refer to section BEFORE YOUR FIRST TREATMENT.

### PULSING AND GLIDING

**Pulsing** is the best technique for working at higher energy levels.

Press the treatment window against your skin. Work in rows (either across or down as shown in Fig 4). Start at one end of each row and progress to the other end. The treatment window creates temporary pressure marks on the treated area so you can see areas you have treated. After each pulse, move the treatment window to the next spot. Keep doing this pulsing action for the area you are treating. **Avoid gaps or overlaps as you use the device.**

**Gliding** is the best technique for working at lower energy levels.

Press the treatment window firmly against your skin and slowly glide over your skin (either across or down as shown in Fig 4). Make sure that the treatment window is always touching your skin. At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes. **Glide at a speed so that you avoid gaps or overlaps.**

- **Do not treat the same area of skin more than once per hair reduction session.**
- **Stop use immediately if you have any discomfort.**

**THE SKIN TONE TESTER**

Light-based hair reduction on darker skin tones may result in adverse effects such as burns, blisters and skin colour changes (hyper- or hypo-pigmentation). The built-in skin tone tester detects the colour of your skin automatically. The skin tone tester checks your skin colour at the start of, and occasionally during, each treatment.

**If your skin tone is too dark for safe use, the device will not operate, and the skin tone tester light will blink.**

**WHAT TO EXPECT DURING USE**

For many people, using the device may be their first experience with a light-based device. During a treatment session it is normal to experience or feel:

- **Noise:** This is the cooling fan
- **Popping sounds:** These occur with each flash
- **A sensation of warmth and tingling:** During each light pulse.
- **Some redness:** It is not uncommon to see some very mild, pink-like skin colouration during or just after a treatment. This is usually most noticeable around the hairs themselves.
- **Stop use immediately if you see your skin is very red.**
- **Contact your physician immediately, if you experience any complication.**

---

**USING YOUR DEVICE**

---

**BEFORE YOUR FIRST TREATMENT**

- Remove all jewellery before use.

You must decide which energy level is right for your skin and the area you want to treat. The device has 5 levels of energy from 1 (the lowest) to 5 (the highest). For example, you might treat one area at energy level 1, and another at level 4. For information on energy levels, refer to section ENERGY LEVELS. To choose the right energy level, test each treatment area separately (refer to the skin colour chart at the front of this manual).

1. Press the control button to turn on the device. The level 1 energy level indicator light will come on.
  2. Place the treatment window on your skin to deliver one pulse. If you do not experience abnormal discomfort press the control button again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.
  3. Continue testing at increased energy levels. If you do not experience any adverse effects within one hour, you can start the full treatment at that energy level.
- **Use a lower energy level if you experience any discomfort.**

Repeat this test for each part of the body you wish to treat.

## YOU ARE NOW READY TO USE YOUR DEVICE

1. Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. **Do not wax, pluck or tweeze the hair out.**
  2. Press the control button. The device will turn on, the level 1 energy indicator light will be on, and a fan sound will start.
  3. Approximately 1 second after pressing the control button, the ready light will turn on. The device is then ready to use at the lowest energy level.
  4. Set the energy level by pressing the control button repeatedly to increase the energy level until you have the level you want. The energy level indicator lights indicate the chosen setting (refer to section ENERGY LEVELS).
  5. Decide if you wish to treat the area by **Gliding** or **Pulsing** (refer to section PULSING AND GLIDING).
  6. Apply the treatment window to the skin. Make sure that the skin is spread evenly and smoothly. When the treatment window is in full contact with the skin, the ready light will start blinking slowly.
  7. The device will first determine the colour of your skin automatically and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.
- **If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the ready light will turn on again, and start blinking slowly when pressed against the skin.**
  - If you are pulsing, repeat the process. If you are gliding, simply continue. Continue to pulse or glide until you have treated the area fully.
  - **Stop use immediately, if your skin blisters or burns.**

## AFTER TREATMENT

When you have finished your treatment session:

1. Turn the device off by giving the control button an extended press. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Remove the adapter from the electrical outlet.

## CLEANING AND MAINTENANCE

It is important to clean your device after each treatment.

1. Unplug the device.
  2. To clean the device including the treatment window, wipe with a clean, dry cloth. Never use acetone or other abrasive fluids, because they will damage your device.
  3. After cleaning, we recommend you store the device in its original box.
- Never immerse your device or any of its parts in water.

## TROUBLESHOOTING

Problem	Check
My device does not start.	Make sure the adapter is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the ready light is on. Make sure that the treatment window is evenly and firmly pressed to the skin and that you are moving the device across your skin. It is also possible that your skin colour is too dark.
The skin tone tester light and the ready light blink.	The skin colour is too dark.
One or more lights blink.	The device has malfunctioned.
The device has stopped working and is in standby mode.	Press the control button again to start the device.

## RECYCLING

To avoid environmental and health problems due to hazardous substances, appliances and rechargeable and non-rechargeable batteries marked with one of these symbols must not be disposed of with unsorted municipal waste. Always dispose of electrical and electronic products and, where applicable, rechargeable and non-rechargeable batteries, at an appropriate official recycling/collection point.

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© 2017 Spectrum Brands, Inc.
Remington® is a Registered Trade Mark of Spectrum Brands,Inc.,
VARTA Consumer Batteries GmbH & Co. KGaA, Alfred-Krupp-Str. 9, 73479 Ellwangen, Germany